

LabPad® INR

My Handy Lab



LabPad® INR



LabPad® INR is an innovative in vitro diagnostics mobile device for PT/INR testing¹. Ideal to monitor patients on Warfarin treatment, LabPad® INR is a connected device designed for strengthening the link with the medical and nursing team.

Thanks to its patented lensless imaging technology² associated with a unique Tsmart® microfluidic microcuvette, LabPad® INR is portable, fast and accurate³. Requiring a very small amount of capillary blood, it is also very easy to use, especially by elderly people⁴.



Features

- User-friendly
- Quick & reliable
- Connected with Bluetooth

A link between patients on the move and healthcare professionals



Intuitive

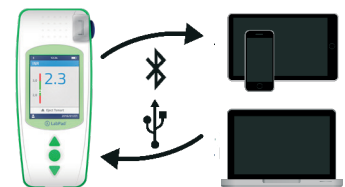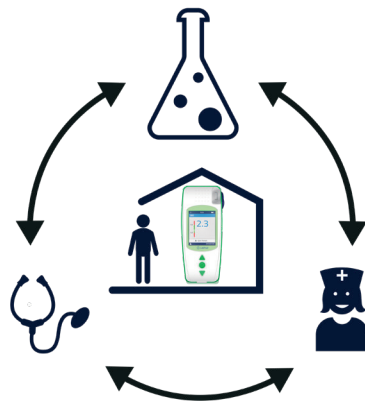

Connected

¹ Garcia-Alamino JM et al, "Self-monitoring and self-management of oral anticoagulation", The Cochrane Library, Issue 5, (2012), p1-66

² Faivre et al., "Coagulation dynamics of a blood sample by multiple scattering analysis", Journal of Biomedical Optics, Volume 16, Issue 5, (2011)

³ Results obtained from a clinical study performed in Grenoble University Hospital on 200 patients (2018); more results obtained through pilot studies are available on demand

⁴ Results obtained from a use study of measurement INR device (150 patients on Warfarin treatment and 100 healthcare professionals), (2012)



A unique T-shaped microcuvette


The Tsmart® microcuvette is specifically designed to ease the drop of capillary blood inside the 3µL central micro channel.

Benefits


- Built-in QC: quick-reading of Datamatrix initiates a series of checks
- Convenient handling: one-way insertion and touchless release
- Storable at room temperature: packed in individual pouch with desiccant

Three simple steps to perform a new test


1

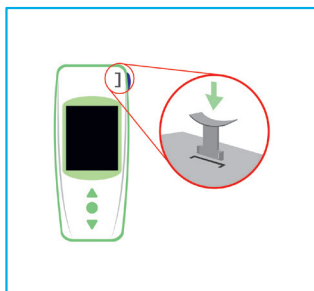

2

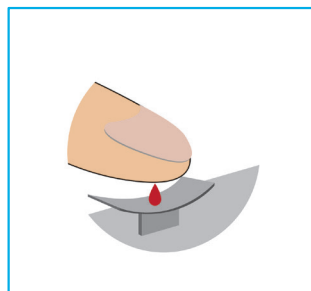

3

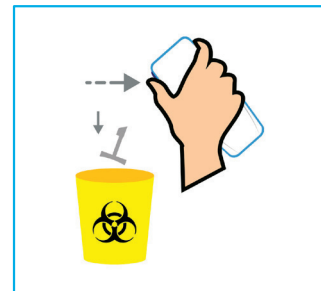

Product Features

3 steps to perform a test
3µL of blood required
3 units displayed: INR, PT and QT
3 buttons to operate

External QC available in 2 levels
Colour screen
Memory capacity of 1,000 results
Built-in battery

Ordering information

LabPad® INR	PT/INR testing device with battery charger, carrying bag and user manuals	03760097080008
Tsmart® INR	box of 12 microcuvettes for PT/INR	03760097080121
	box of 24 microcuvettes for PT/INR	03760097080138
	box of 48 microcuvettes for PT/INR	03760097080145

Avalun SAS
7, Parvis Louis Néel
38000 Grenoble
FRANCE

©2018

For further information visit the website www.mylabpad.com
Photos credit : Avalun company

