HemoSmart [™]*GOLD* ■ Hemoglobin Screening System 13:30 03/23 HemoSmart GOLD **Precious metal** electrode for optimum accuracy ApexBio Hgb/Hct Hgb / Hct display Quick result Splash proof Strip ejection 1µL **APEXBIO**


HemoSmart GOLD Hemoglobin Screening System

HemoSmart GOLD offers a fast and easy way to measure hemoglobin level and is the first biosensor-based palm sized hemoglobin testing system in the world. The one of its kind strip-based testing method offers ease of use and portability and requires only 1µL of blood. Testing is done away from the instrument so cross contamination can be minimized. The instrument is made of splash proof material and can withstand rugged cleaning by chlorine-base disinfectants. The large, back-lit LCD provides big and clear readout even in poor lighting conditions. Results can be easily upload into laptop computers or worksation via USB connection.

Easy to use





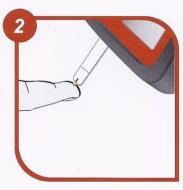

Apply sample



Results in 5 seconds

Specifications

■ Measurement Range: 7.0-24.0 g/dL

Quick Reaction Time: 5 seconds

■ Tiny Blood Sample: 1 µL

■ Operating Temperature: 10 °C~40 °C


Memory: 100 sets with date and time

Power Supply: 3A battery X 2

Additional Features:

LCD backlight **USB** connection Bluetooth (Optional)

Accuracy and Precision Study vs. Reference

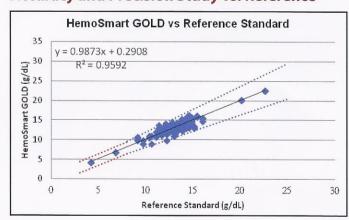

Number of reading	100	100	100	100	100
Mean (g/dL):	7.6	10.3	12.9	16.6	21.8
S.D. (g/dL):	0.3	0.4	0.4	0.5	0.3
C.V. %:	3.5	3.8	2.9	2.7	1.5

No.7, Li-Hsin Rd. V, Hsinchu Science Park, Hsinchu, Taiwan, R.O.C.

TEL: 886-3-5641952 / FAX: 866-3-5678302 / E-Mail: info@apexbio.com / Website: www.apexbio.com