




COMPRESSOR NEBULIZER NE-C25S

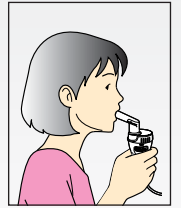

Using the Mouthpiece

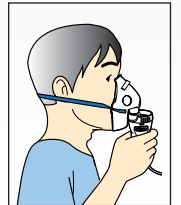

Using the Adult Mask

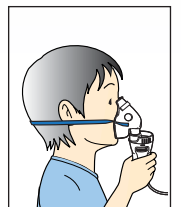

Using the Child Mask

FOR EFFICIENT RESPIRATORY THERAPY

- Easy to operate
- High nebulization rate ensures optimized treatment duration
- Efficient delivery of medication
- Compact & travel convenient


Compressor Nebulizer NE-C25S

OMRON Compressor Nebulizer is developed for efficient treatment of asthma, allergies and other respiratory disorders.

EFFICIENT NEBULIZER KIT

The nebulizer kit provides a dense aerosol in the therapeutic range. It has a 6ml medicine cup and leaves minimal residual medication.

The compressed air from the compressor is forced through the nozzle, as it does so, it is mixed with the medication in the baffle and forced upwards. The medication is then turned into a fine aerosol and directed upwards by the baffle, where it is broken down into many smaller aerosol particles.

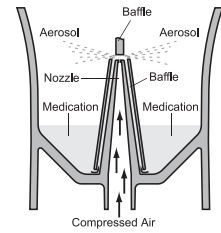







SPECIFICATIONS

Model	NE-C25S
Type	Compressor Nebulizer
Rating	AC 230, 50Hz 1A
Power Consumption	75W
Nebulization Rate	Approx 0.25ml/min (by weight loss)
Particle Size MMAD	4.21 um MMAD = Mass Median Aerodynamic Diameter
Medication Tank Capacity	6ml maximum
Appropriate Medication Quantities	2ml minimum – 6ml maximum
Operating Temperature / Humidity	+10°C to +40.0°C / 30% to 85% RH
Storage and Transport Temperature / Humidity / Air Pressure	-20°C to +60°C / 10% to 95% RH / 700-1060hPA
Weight	Approx. 1.5kg (compressor only)
Dimensions	Approx. 158 (W) × 175 (D) × 100 (H) mm (compressor only)
Contents	Compressor, Nebulizer Kit, Air Tube (180cm), Mouthpiece, Child Mask, 5 pcs Spare Air Filters, Carrying Case, Instruction Manual.

NOTE: These specifications are subject to change without notice.

This OMRON product is produced under the strict quality system of OMRON HEALTHCARE Co., Ltd., Japan.

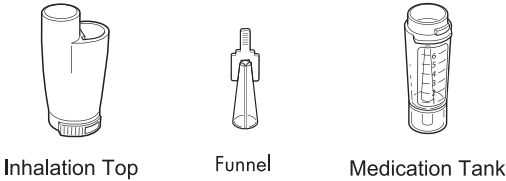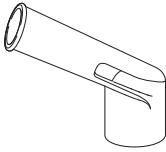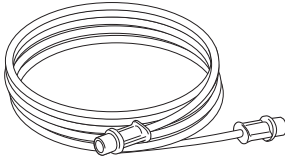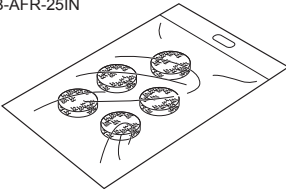Accessories (Sold seperately)

Nebulizer Kit x 1 NEB-NK-25IN 		
Mouthpiece x 1 NEB-MP-25IN 	Air Tube (180cm) X1 NEB-TP-25IN 	
Spare Air Filters x 5 NEB-AFR-25IN 	Child Mask X 1 NEB-MSMP-25IN 	Adult Mask X 1 NEB-MSLP-25IN 



OMRON HEALTHCARE INDIA PVT. LTD.

Office No. 103-106 & 108, 1st Floor,
 Sewa Corporate Park, MG Road, Gurgaon, Haryana - 122001, India
 Email : callcentre@ap.omron.com • Website : www.omronhealthcare.in
Toll Free : 1800-419-0492 - Monday to Friday - 9 am to 6 pm (Except Holidays)

Distributed by: