



OPTIMA POCKET DIAGNOSTIC SETS

TIMESCO OPTIMA POCKET DIAGNOSTIC SETS



The Timesco Optima Pocket diagnostic sets have been designed to offer the Clinician and student alike with a perfect diagnostic tool for ophthalmology and aural examination with the convenience of being small enough to be carried in a pocket!

The Timesco Optima Pocket diagnostic sets feature the same superb quality optics, precision lenses, fibre optics and Xenon illumination for pure white light and durable materials as our desk and wall mounted Diagnostic sets.

Individual Ophthalmoscopes and Otoscopes as well as a combination of both in Diagnostic sets in hard and soft cases are available.

The complete Optima Pocket Diagnostic sets are guaranteed for a period of two years from date of purchase for faults in manufacture and materials. Bulbs, speculums, and any consumables are not covered by guarantee.



OPTIMA POCKET OPHTHALMOSCOPE

The Optima Pocket Ophthalmoscopes feature durable construction, superb bright Xenon white light illumination and five apertures: large, small macular (spot), half moon, fixation and red free filter.

Apertures:



Practitioners can determine the clinical application for each aperture.

Standard spot: for dilated pupil, general application.

Small spot: for undiluted pupils, ideal aperture to start visualization.

Half-moon: for small reflection with small pupils.

Fixation aperture: for observation of eccentric fixation.

Red free aperture: for observation of defined retinal abnormalities.

Lenses:

There are 18 colour coded lenses, green for + and red for -.

+	1	2	3	4	6	8	10	15	20
-	1	2	3	4	6	8	10	15	20

Optima Pocket diagnostic sets in soft pouches are available with and without spare bulbs to make these even less expensive.



2505.250.25
Optima Pocket Ophthalmoscope
In Pouch With 2.5V Spare Bulb

OPTIMA POCKET OTOSCOPE

The Optima Pocket Oscopes features high intensity Xenon illumination and fiber optics and are constructed from durable plastics, metal alloys and stainless steels.

The Oscope can be used with the Optima disposable and reusable speculums which are available in various diameters.

A 3X magnification swivel lens is incorporated in the otoscope proximal head which when rotated allows instruments to be used for aural examinations. Without the speculums mounted the otoscope head can also be used as an oral illuminator.

An insufflation port is also provided in the Oscope head for connection to insufflation bellows for pneumatic testing of the tympanic membrane.

Optima Pocket diagnostic sets in soft pouches are available with and without spare bulbs to make these even less expensive.



2505.100.25
Optima Pocket Oscope In Pouch With
2.5V Spare Bulb

OPTIMA POCKET SETS IN HARD CASE



2505.275.75
Optima Pocket Diagnostic Set In Hard Case With
2.5V Spare Bulbs



2505.250.50
Optima Pocket Ophthalmoscope In Hard
Case with 2.5V Spare Bulb




2505.260.50
Optima Pocket Oscope In Hard
Case with 2.5V Spare Bulb

OPTIMA POCKET SET ACCESSORIES



Optima Pocket Otoscope Speculums Disposable (2.5mm and 4.0mm)

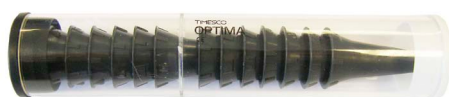

2505.900.25 Optima Pocket Otoscope Speculums Disposable In Tube, 6 Each 2.5mm and 4mm



2505.950.50 Optima Insufflation Tube



2505.325.10 Optima Pocket Otoscope Bulb Only, 2.5V Pack of 6



2505.325.20 Optima Pocket Ophthalmoscope Bulb Only, 2.5V Pack of 6



2505.950.25 Optima Otoscope Insufflation Connector

OPTIMA POCKET DIAGNOSTIC PRODUCT LIST

CAT. NO.	DESCRIPTION
2505.100.25	Optima Pocket Otoscope In Pouch With 2.5V Spare Bulb
2505.105.25	Optima Pocket Otoscope In Pouch Without 2.5V Spare Bulb
2505.250.25	Optima Pocket Ophthalmoscope In Pouch With 2.5V Spare Bulb
2505.255.25	Optima Pocket Ophthalmoscope In Pouch Without 2.5V Spare Bulb
2505.250.50	Optima Pocket Ophthalmoscope In Hard Case With 2.5V Spare Bulb
2505.260.50	Optima Pocket Otoscope In Hard Case With 2.5V Spare Bulb
2505.275.75	Optima Pocket Diagnostic Set In Hard Case With 2.5V Spare Bulbs
2505.300.10	Optima Pocket Ophthalmoscope Head And 2.5V Bulb Only
2505.300.20	Optima Pocket Otoscope Head And 2.5V Bulb Only
2505.300.30	Optima Pocket Handle Only
2505.325.10	Optima Pocket Otoscope Bulb Only, 2.5V, Pkt 6
2505.325.20	Optima Pocket Ophthalmoscope Bulb Only, 2.5V, Pkt 6
2505.850.25	Optima Pocket Otoscope Spec Disposable 2.5mm, Pack of 250
2505.850.40	Optima Pocket Otoscope Spec Disposable 4.0mm, Pack of 250
2505.900.25	Optima Pocket Otoscope Spec Disposable In Tube, 6 Each 2.5mm and 4mm
2505.905.25	Optima Pocket Soft Pouch For Otoscope And Ophthalmoscope
2505.910.25	Optima Pocket Hard Case For Diagnostic Set
2505.915.25	Optima Pocket Hard Case For Ophthalmoscope
2505.920.25	Optima Pocket Hard Case For Otoscope
2505.950.25	Optima Otoscope Insufflation Connector
2505.950.50	Optima Insufflation Tube



Timesco Healthcare Ltd, Timesco House, 3 Carnival Park, Carnival Close, Basildon, Essex, SS14 3WN, England

T +44 (0) 1268 297 700 | F +44 (0) 1268 297 800

E info@timesco.com | www.timesco.com

Ref: B.D.010/2. E & OE April 2014.

TIMESCO
INSTRUMENTS FOR LIFE



TIMESCO ENGLAND