



**üzümcü**



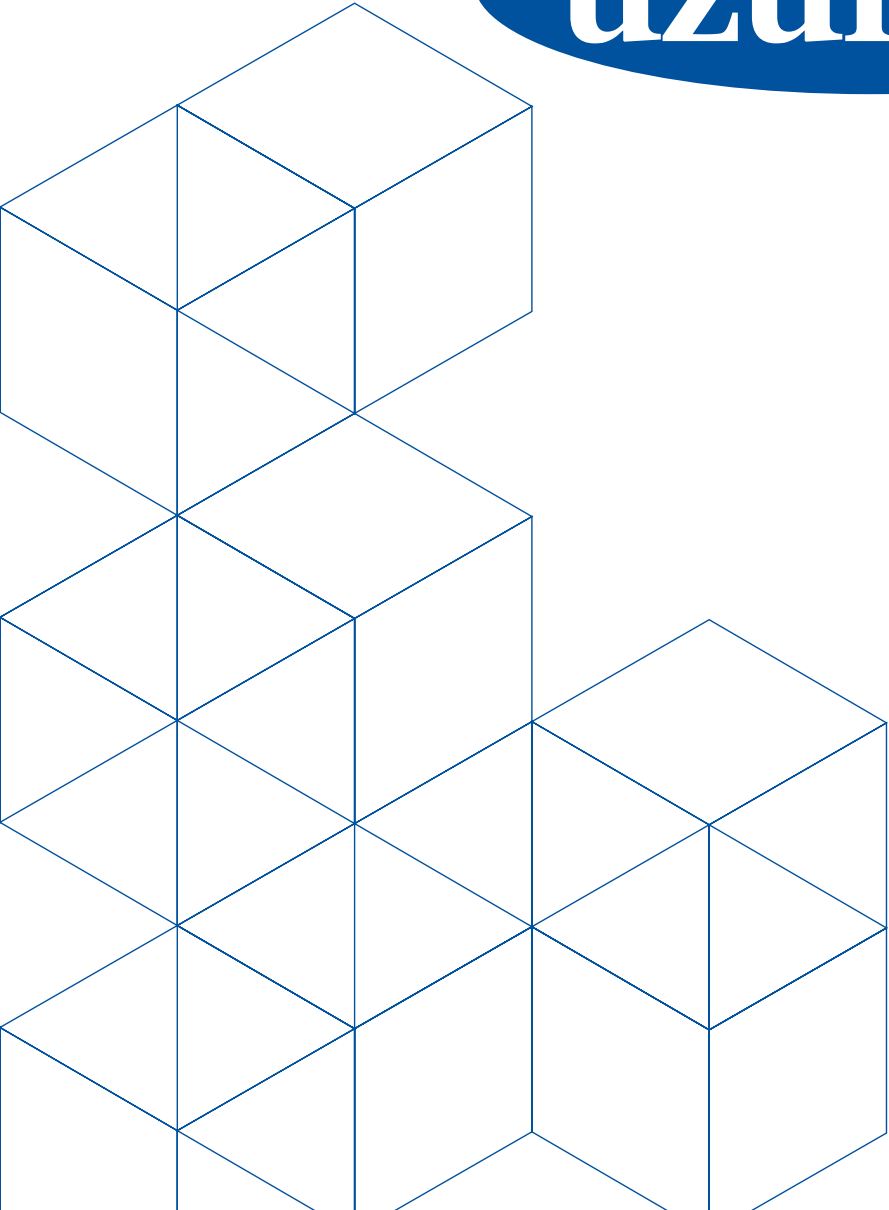

üzümcü

# OPERATING TABLES & ACCESSORIES

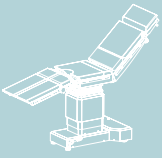

OPERATING  
TABLES  
&  
ACCESSORIES



We are better together  
We enable you to treat your patients  
on the best tables.

üzümcü

# OPERATING TABLES

## Electrical Movements

- Height
- Back Rest
- Trendelenburg
- Reverse Trendelenburg
- Lateral Tilt
- Sliding

## Mechanical Movements


- Head Rest
- Kidney Surgery Position
- Leg Rests, Up / Down
- Leg Rests, Right / Left
- Central Locking

## Sizes

- Width x Length : 610x2130 mm (2P)  
: 610x2040 mm (2V)
  - Height (min) : 810 mm (2P) | 840 mm (2V)
  - Height (max) : 1210 mm (2P) | 1240 mm (2V)
- \*Max. & min. dimensions are measured with mattress.*

Üzümcü OP-2P operating tables are designs which can be adjusted to get positions suitable for carrying out all surgical operations, produced by considering the patient and user requirements and maximize the operational efficiency.



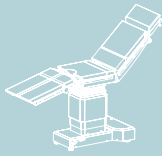

## Angles

• Sliding Movement	: 300 mm	• Leg Rest Up / Down	: +65°/-90°
• Lateral Tilt Movement	: +20°/- 20°	• Leg Rest Parallel Movement	: 80°
• Back Rest Movement	: +90°/- 40°	• Head Rest Movement	: +30° / - 60°
• Trendelenburg	: +25°	• Kidney Surgery Position	: 250°
• Reverse Trendelenburg	: -25°	• Flex / Reflex Position	: +220° / -90°

The optional patient heating system of the OP-2V table enables to keep the patient's heat stable with 0.1°C accuracy. The system shows measured and set temperatures on a LCD display at the same time. The heating device has an advanced alarm and safety system.

• Three layers visco mattress for ensuring patients comfort

• Positioning with Hand, Body, and Remote Control.

• Telescopic Height Adjustment



• Patient heating system

< OP-2V



# OPERATING TABLES

## Electrical Movements

- Height
- Back Rest
- Trendelenburg
- Reverse Trendelenburg
- Lateral Tilt
- Sliding

## Mechanical Movements


- Head Rest
- Kidney Surgery Position
- Leg Rests, Up / Down
- Leg Rests, Right / Left
- Central Locking

## Sizes

- Width x Length : 610x2200 mm (6N)  
: 780x2140 mm (6L)
  - Height (min) : 800 mm (6N) | 820 mm (6L)
  - Height (max) : 1200 mm (6N) | 1220 mm (6L)
- \*Max. & min. dimensions are measured with mattress.*



# OM-6N / OM-6L

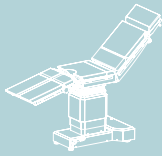

OPERATING  
TABLES  
& ACCESSORIES

## Angles

• Sliding Movement	: 300 mm (6N)	• Leg Rest Up / Down	: +65°/-85°
• Lateral Tilt Movement	: +20°/- 20° (6N)   +30°/- 30° (6L)	• Leg Rest Parallel Movement	: 90° (6N)   80° (6L)
• Back Rest Movement	: +90°/- 40°	• Head Rest Movement	: +35° / - 80° (6N)   +35°/-60° (6L)
• Trendelenburg	: +25°	• Kidney Surgery Position	: 250°
• Reverse Trendelenburg	: -25°	• Flex / Reflex Position	: 210° / 115° (6N)   220°/90° (6L)



Nano surface protection technology is easy to clean, highly resistant to corrosion and has an antibacterial effect that prevents the production and inhabitation of blood and body fluids on the surface.



# OPERATING TABLES

## Detailed Specifications



### < Functional Design & Wide Range of Accessory Options

- Modular design taking shape by different surgical operations with high stability
- Solutions ensuring to reach the position required by any kind of surgical operations with its wide range of accessory options
- Telescopic body specifications enabling height adjustment at different levels
- Wheel system which is made of ABS plastic, with high multi-directional mobility and 125 mm in diameter (OP Series)



### < Sliding

- Image reception ability in a wide area without disturbing the sterile environment with the flouroscopy device by its slideability within 300 mm on horizontal plane



### < Optimal Clearance In Surgical Area

- Maximum surgical area visibility and operational success with kidney surgical elevator integrated on the table and that can be adjusted up to 300°(120 mm)



### < Different Mattress Options


- High-quality image reception facility with table top made of compact material with high x-ray permeability
- Visco or polyurethane mattress options ensuring patient comfort by taking shape of body parts of patient and preventing decubitus prophylaxis
- Cleanable and sterilizable mattress for decontamination of operating room with its seamless structure and demountability



### < Modular Design

- Design that can get different angles according to the operation's kind, demountable head and foot plates
- Unlimited perineal accessibility with demountable foot plate



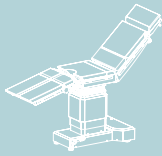

## < Easy Approach For C-Arm

- Possibility to have high image quality from all parts of the table with X-ray permeable table top
- Suitable base structure for easy access and approach of the C-arm scopy and surgical microscopes



## < Accessory Positioning and Battery System

- Attachment system mountable on stainless steel rail integrated on both sides of the table, enabling connecting a new accessory without demounting existing accessories and preventing loss of time.
- Possibility of positioning the table in case of electricity cuts thanks to the rechargeable battery system that guarantees thirty movements.




## < T Base Specifications

- T Base enabling surgical team and equipment to approach the table optimally and contributing to operational success
- 316 quality stainless steel body resistant to fractures and corrosion
- Silent and nonvibrating operation



## < Remote Control

- Multi control system with wireless touch screen tablet, hand set and body control panel
- Rapid recall in case of need by storing five different surgery positions
- Position reset by a single button
- Easy use in endosurgery operations thanks to the illuminated screen

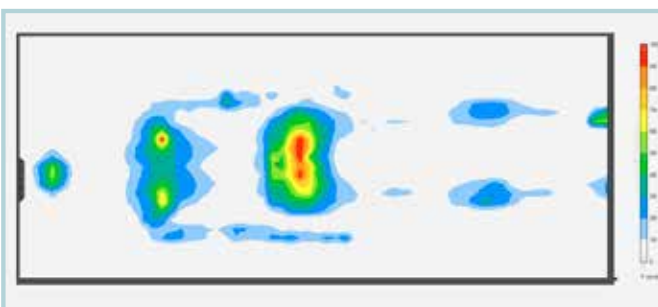

## < Optimum Mattress Performance

- Uzumcu operating table mattresses are designed and manufactured according to the Pressure Mapping test results. These tests are performed in house test systems in Uzumcu factory.

# OPERATING TABLES

## Mechanical Movements

- Height Adjustment
- Back Rest
- Trendelenburg
- Reverse Trendelenburg
- Lateral Tilt Movement
- Head Rest Movement
- Leg Rests, Up / Down
- Leg Rests, Right / Left

## Sizes

- Width x Length : 610 x 2150 mm
  - Height (min) : 735 mm
  - Height (max) : 985 mm
- \*Max. & min. dimensions are measured with mattress.*




Üzümçü OM-2M table offers suitable solutions for every surgery branch and budget, it is a model easy to use, long-lasting with low cost and enables all positions mechanically.



*\*Electrical height adjustment in telescopic body (OM-2ME)*



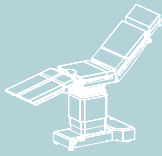

## Angles

• Lateral Tilt Movement	: +30°/- 30°	• Leg Rest Up / Down	: +60°/-85°
• Back Rest Movement	: +70°/- 35°	• Leg Rest Parallel Movement	: 90°
• Trendelenburg	: +20°	• Head Rest Movement	: +35° / -80°
• Reverse Trendelenburg	: -25°	• Flex / Reflex	: 215°/110°

## Detailed Specifications



### < Functional Design & Wide Range of Accessory Options

- Modular design taking shape by different surgery operations with high stability
- Solutions ensuring to reach the position required by any kind of surgery operations with its wide range of accessory options
- Telescopic body enabling height adjustment at different levels



### < Ergonomic Positioning Arms

- Trendelenburg/reverse trendelenburg and spine positioning arms that are easily accessible and enables ergonomic use



### < Table Top Specifications and Lateral Use

- Cleanability and sterilizability suitable for decontamination with its seamless structure and demountability
- High-quality image reception facility, table top being made of compact material with high x-ray permeability
- Right-left lateral tilt arm, practical to use and integrated with the operation table, enables mobility for the surgery team and can be hidden under table sitting plate when not in use



### < T-Base Specifications

- T-Base contributes to operational success by enabling optimal approach of surgery team and equipment to table
- 316 quality stainless steel body resistant to fractures and corrosion
- Silent and non-vibrating operation



### < Breaking System & Hydraulic Lever

- Wheel system enabling mobility in the operation room
- Brake system enabling positioning in the required position
- Telescopic body with height adjustment at required level by hydraulic pedal

# OPERATING TABLES



Trendelenburg Position



Reverse Trendelenburg Position



Left Lateral Tilt Position



Right Lateral Tilt Position



Reflex Position



Flex Position



Chest (Breast) Surgery Position



Struma Surgery Position



Bariatric Laparoscopic Surgery Position



Anaesthesia Preparation and General Surgery Position




Chest Surgery Position



Kidney, Sternum Surgery Position

# POSITIONS

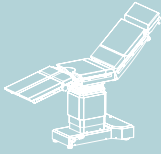

OPERATING  
TABLES  
& ACCESSORIES



Kidney Surgery Position



Gynecologic and Urologic  
Surgery Position



Gynecologic, Urologic,  
Colorectal Surgery Position



Spinal, Proctologic Surgery  
Position



Spinal Surgery Position



Neurosurgery Position



Neurosurgery Position



Neurosurgery Position



Arthroscopic Knee Surgery Position



Meniscus Surgery Position



Orthopedic Fracture Surgery  
Position



Ophthalmology Position

# OPERATING TABLES



OM-120 Arm Support  
AC: 11.1020



OM-110 IV Stand  
OM-100 Anaesthesia Screen\*  
AC: 11.1019 \*AC: 11.1018



OM-150 Body Strap  
AC: 11.1021



OM-220 Geopel Knee Crutch (Pair)  
AC: 11.1030



OM-600 Leg Positioning Boot  
AC: 11.1063



OM-380 Irrigation Bowl  
AC: 11.1032



OM-152 Hand-Wrist Strap  
AC: 11.1023




OM-155 Foot-Wrist Strap  
AC: 11.1024



OM-405 Clamp  
AC: 11.1105

# ACCESSORIES

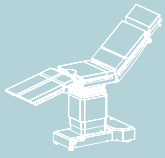

OPERATING  
TABLES  
& ACCESSORIES



OM-390 Clamp  
A.C.: 11.1027



OM-130 Shoulder Support (Pair)  
A.C.: 11.1022




OM-190 Lateral Support (Pair)  
A.C.: 11.1025



OM-200 Foot Support (Pair)  
A.C.: 11.1026



OM-310 Proctological Accessories  
A.C.: 11.1033

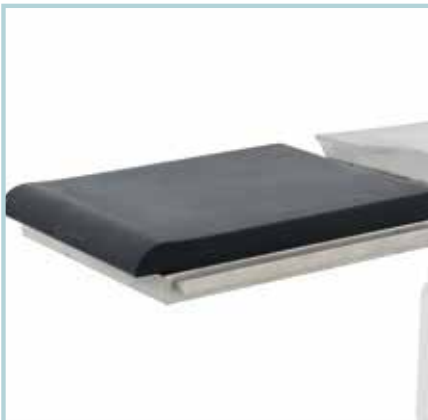

OM-255 Extension Leg Plate (Single)  
A.C.: 11.1039



OM-330 X-Ray Cassette Holder  
A.C.: 11.1040



OM-990 Accessory Trolley  
A.C.: 11.1043

# OPERATING TABLES



OM-240 Knee Arthroscopy Support  
A.C.: 11.1037



OM-250 Meniscus Positioning Roll  
A.C.: 11.1038



OM-520 Orthopedic Traction Set  
A.C.: 11.1042



OM-365 Arm Surgery Plate  
A.C.: 10.1041



OM-550 Orthopedic Traction Set  
Accessory Trolley  
A.C.: 11.1044



OM-270 Ophthalmologic Head Support  
A.C.: 11.1036



OM-285 Neurosurgery Head Rest With  
Pin & Horseshoe Head Rest  
A.C.: 11.1034




OM-300 Croos Bar  
A.C.: 11.1035



OM-460 Mattress, Adult, Gel  
A.C.: 11.1050  
OM-450 Mattress, Pediatric, Gel




# ACCESSORIES

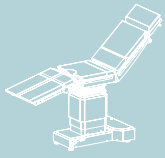


OPERATING  
TABLES  
& ACCESSORIES




OM-428 Foot Wrist Rest, Gel  
A.C.: 11.1011

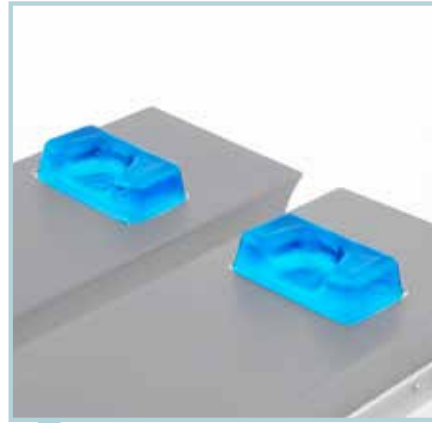


OM-425 Heel Rest, Gel  
A.C.: 11.1010

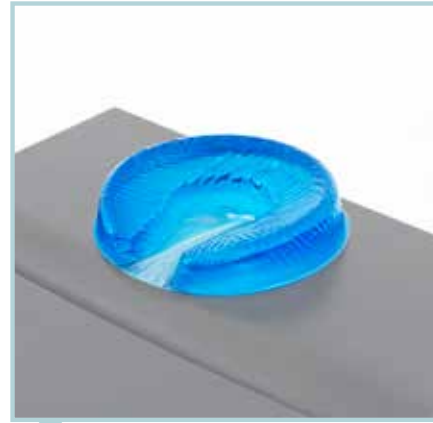

OM-410 Horseshoe Head Rest, Gel  
OM-412 Horseshoe Head Rest, Gel, Pediatric\*  
A.C.: 11.1007 \*A.C.: 11.1008

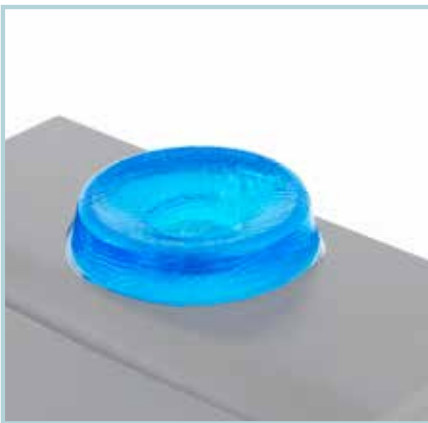


OM-430 Head Rest, Gel  
OM-432 Head Rest, Gel, Pediatric\*  
A.C.: 11.1012 \*A.C.: 11.1013

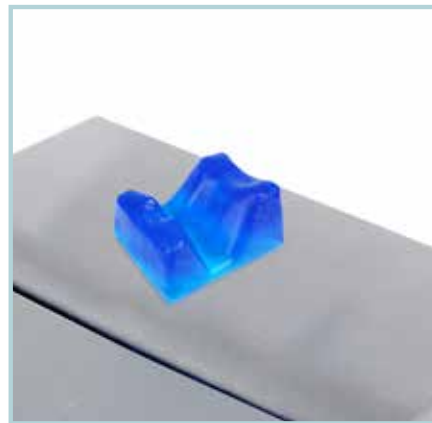

OM-415 Prone Head Support  
A.C.: 11.1102




OM-420 Arm Rest, Gel  
A.C.: 11.1009

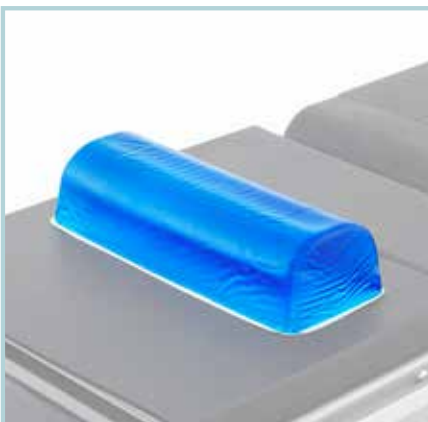

OM-440 Chest Cushion, Gel  
OM-442 Chest Cushion, Geli Pediatric\*  
A.C.: 11.1014 \*A.C.: 11.1015



OM-446 Spinal Surgery cushion with silicone  
OM-447 \*Spinal Surgery cushion without silicone  
A.C.: 11.1103 \*A.C.: 11.1104

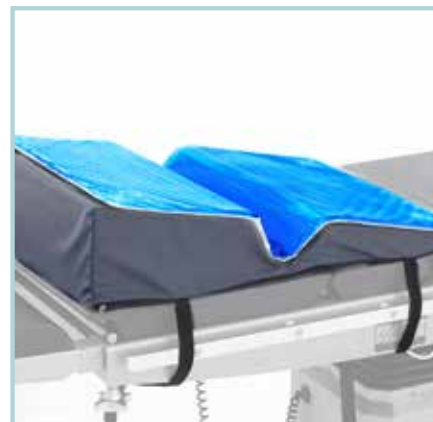

OM-445 Silicon lateral position support  
A.C.: 11.1016

# C-ARM IMAGING TABLE

## Table Movements

- Height
- Multifunctional Table Top Slide


## Mechanical Movements

- Central Locking

## Sizes


- Width x Length : 465x2450 mm
- Height (min) : 845 mm
- Height (max) : 1105 mm

Üzümcü XT-15 imaging tables are models that enable excellent C-arm accessibility with non-metallic carbon fiber top tray and maximum imaging quality at low radiation, offering solutions for surgery and interventional radiology.

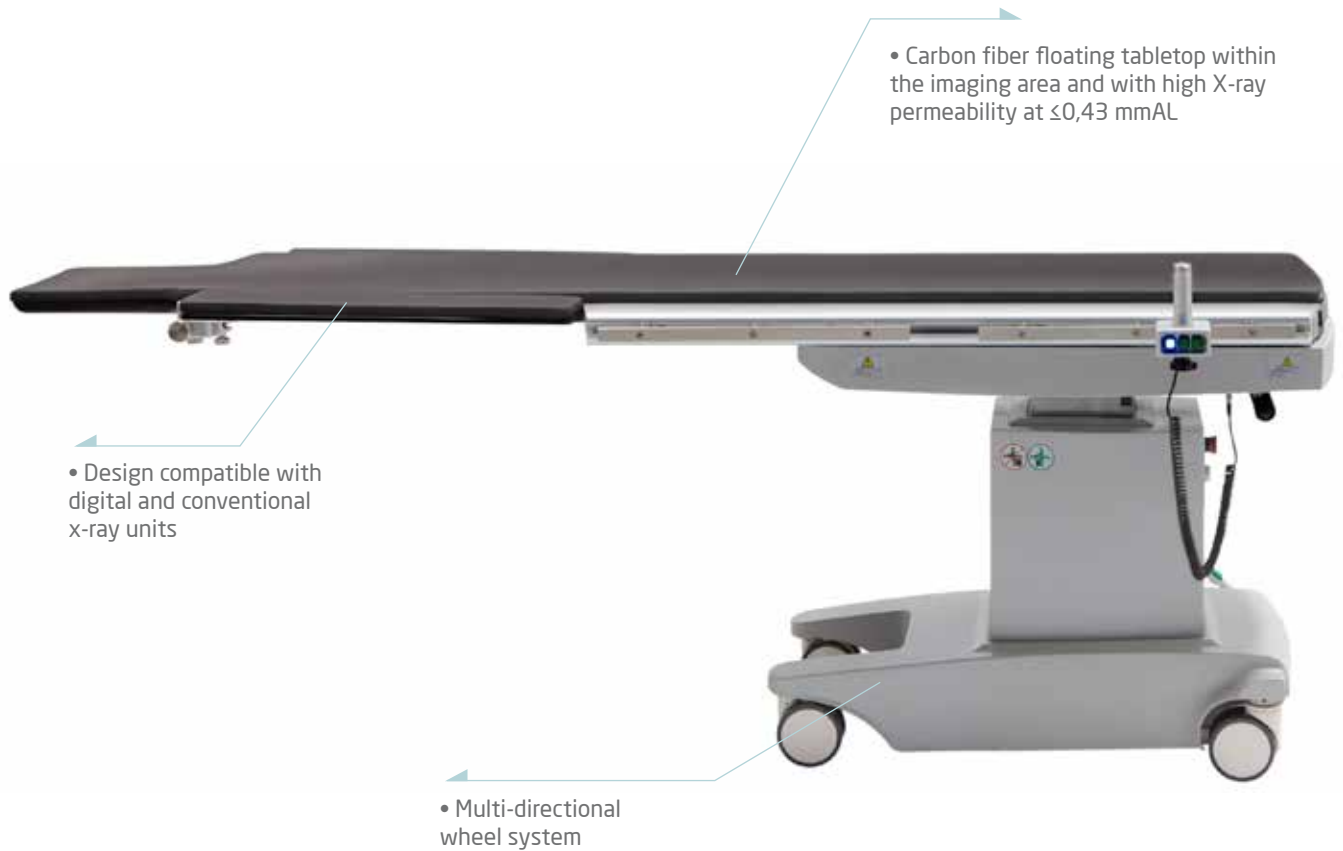

XT-15 ^



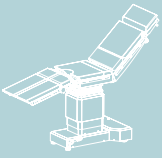

## Table Top Movements

- X Axis Movement : 820 mm
- Y Axis Movement : 220 mm
- Vertical Movement : 275 mm

1100 mm carbonfiber tabletop without any metal part enables various patient positionings. Low X-Ray absorption carbonfiber materials get high clear images. It extends the X-Ray machine Spheroid's tube lifetime.



- Ergonomic hand control unit enables height adjustment and floating plate movements



# GYNECOLOGICAL TABLE

## Table Movements

• Height Adjustment	: 720 - 970 mm	• Trendelenburg	: 20°
• Back Rest	: 50°	<i>*Max. &amp; min. dimensions are measured with mattress.</i>	

• Visco elastic seamless mattress ensuring hygiene conditions and patient comfort. It can be easily cleaned in accordance with hospital decontamination rules

• Easily replaceable leg and foot rest accessories can be used for different gynecological operations

• Durable steel body with electrostatic paint


Üzümcü JM-4F Gynecological Table is a model suitable for any kind of gynecological examination and operation, contributing visually to the operation environment with its different colour options.

• Ergonomic examination and diagnose possibility for the doctor

< JM-4F

• Height adjustment, Trendelenburg position, and back section can be adjusted both by hand and foot control for easy and fast positioning



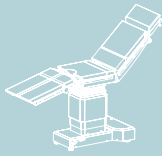

## Detailed Specifications



### < Waste Collection

- Sustainable hygiene conditions with waste collection system, which is easy to clean and disinfect thanks to the round concave structure



### < Hygiene

- Roll cover dispenser at spine plate for patient's comfort and hygiene requirement



### < Comfort

- Table extension which is integrated within the sitting plate frame and moveable by telescopic rail system in case of need

# üzümcü

## ÜZÜMCÜ Tıbbi Cihaz ve Medikal Gaz Sistemleri San. ve Tic. A.Ş.

### Headquarters & Factory

Ankara Konya Devlet Karayolu, 29. Km.

06830 Gölbaşı / Ankara TURKEY

p: +90.312.615 53 53

f: +90.312.615.53 90

### İstanbul Regional Headquarters

Oruç Reis Mahallesi Tekstil Kent Cad.

Tekstil Kent Koza Plaza B Blok - Kat 10

34235 Esenler / İstanbul, TURKEY

p: +90.212.438 39 00

f: +90.212.438 39 49

[uzumcu@uzumcu.com.tr](mailto:uzumcu@uzumcu.com.tr)

[www.uzumcu.com.tr](http://www.uzumcu.com.tr)



KT-TB-20.01.16/3

All information in this document is subject to change without prior notice. No parts of this document may be reproduced or transmitted in any form by any means without the express written permission of the manufacturer. All brand name and product names in this document are trademarks or registered trademarks of ÜZÜMCÜ