

KIT FOR THORACENTESIS AND PARACENTESIS




THORACENTESIS

- A practical system with single-acting valves
- Available in various configurations
- Wide range of needles

PARACENTESIS

- Available with bags in various capacities
- Atraumatic cannula

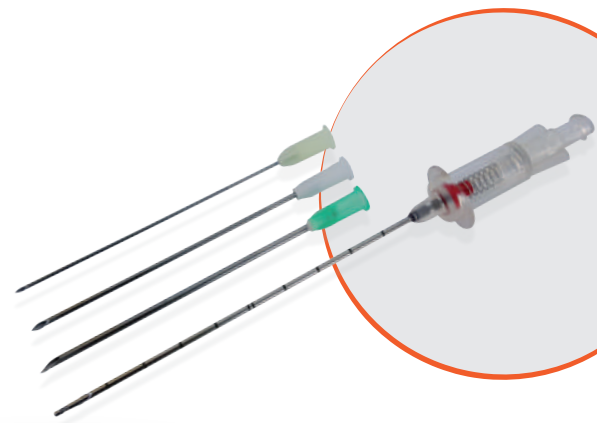

KIT FOR THORACENTESIS

The device is used for aspiration of pleuric exudate.

The kit contains everything necessary for performing thoracentesis (needles, syringe and collection bag).




The aspiration set is highly versatile in use due to availability of different needle configurations which differ in dimension and operating characteristics.

- Trio of 55 mm needles available in the diameters 14G, 16G and 19G
- Trio of 80 mm needles available in the diameters 14G, 16G and 19G
- Veress needle with atraumatic retractable tip.

The "TORAMATIC" device allows safe aspiration as a result of being equipped with one-way valves which prevent re-infusion of any of the aspirated fluid back to the patient.

As an alternative to the "TORAMATIC" Bioservice offers the "TORASET" version, equipped with a tightly sealed 3-way stopcock with directional indicators of the fluid aspiration and subsequent drainage into the bag.

The device is supplied in EtO sterilized packaging

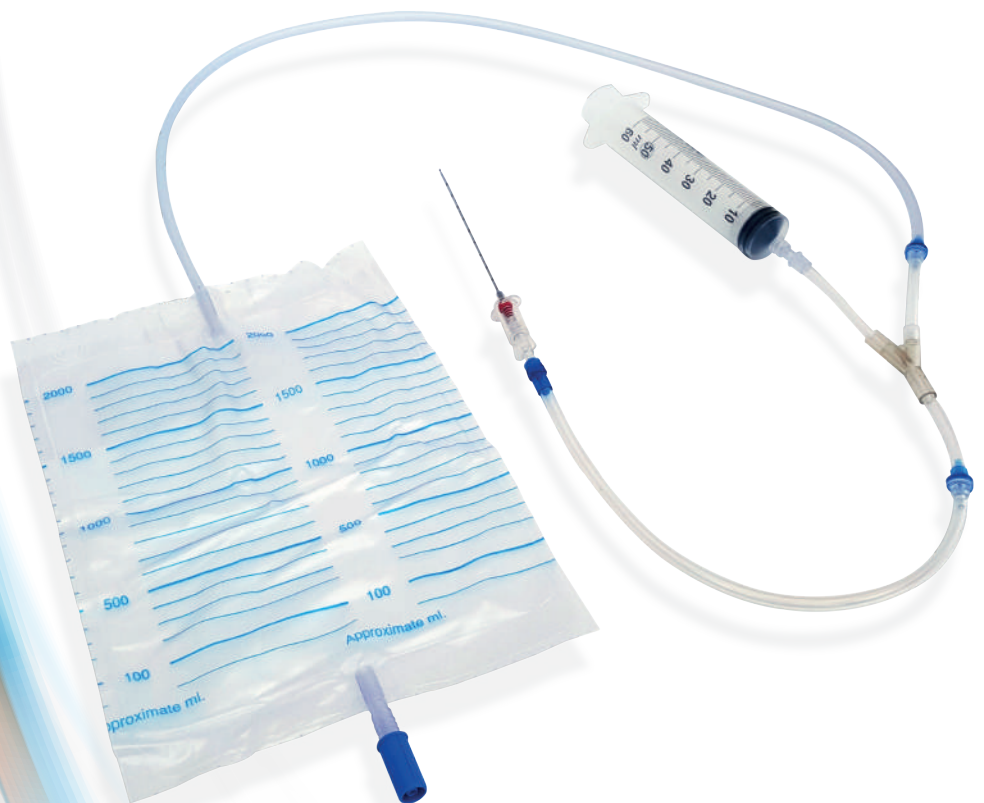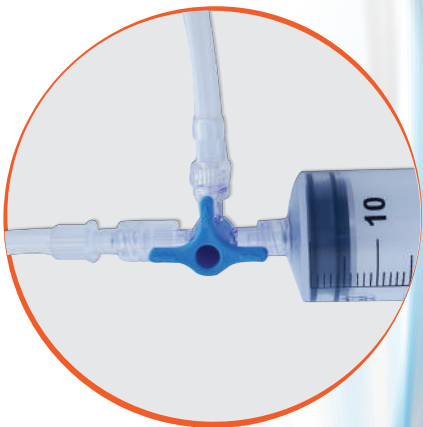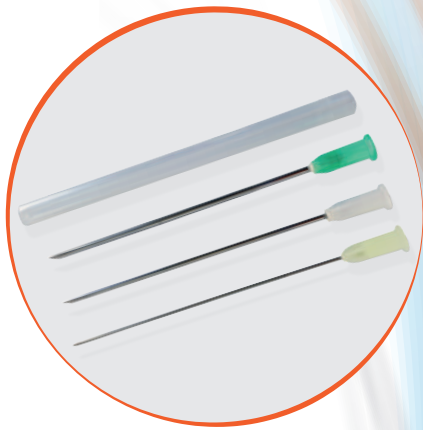

KIT FOR PARACENTESIS

The device is used for aspiration of ascitic fluid.

The kit contains every thing necessary for performing paracentesis (needles, syringe and collection bag).

The set is available with graduated bags in different capacities:

- 2-litre bag (Model PARASET 2)
- 5-litre bag (Model PARASET 5)
- 8-litre bag (Model PARASET 8)
- 12-litre bag (Model PARASET12)

All the bags are fitted with a non-return valve and drainage tap. The device is equipped with a tightly sealed 3-way stopcock with directional indicators of the fluid aspiration and subsequent drainage into the bag.

"PARASET" contains a cannula needle which ensures that the paracentesis procedure is less traumatic for the patient.

After perforation and subsequent extraction of the internal metal cannula, the flexible catheter which remains inserted in the patient prevents any type of internal trauma.

The catheter reassumes its original form after any accidental bending, thus preventing obstruction of the drainage fluid passage.

On request of the customer, cannula needles with different dimensions are available.

The kit can also be made up with a Veress needle.

The device is supplied in EtO sterilized packaging.



KIT FOR THORACENTESIS AND PARACENTESIS



KIT FOR THORACENTESIS

| CODE | MODEL | CHARACTERISTICS | PACKAGING |
|---------|-----------------|---|-----------|
| H9200E | TORASET 1-55 | - trio of 55 mm needles with protective cap - 3-way stopcock - 60 ml syringe - 2-litre drainage bag with non-return valve and drainage tap | 30 pcs |
| H92010E | TORASET 1-80 | - trio of 80 mm needles with protective cap - 3-way stopcock - 60 ml syringe - 2-litre drainage bag with non-return valve and drainage tap | 30 pcs |
| H92020E | TORASET 1-VER | - 100 mm Veress needle - 3-way stopcock - 60 ml syringe - 2-litre drainage bag with non-return valve and drainage tap | 30 pcs |
| H92030E | TORASET 2-55 | - trio of 55 mm needles with protective cap - 3-way stopcock - 60 ml syringe - 2-litre drainage bag | 30 pcs |
| H92040E | TORASET 2-80 | - trio of 80 mm needles with protective cap - 3-way stopcock - 60 ml syringe - 2-litre drainage bag | 30 pcs |
| H92050E | TORASET 2-VER | - 100 mm Veress needle - 3-way stopcock - 60 ml syringe - 2-litre drainage bag | 30 pcs |
| H92001E | TORAMATIC 1-55 | - trio of 55 mm needles with protective cap - 2 one-way valves - 60 ml syringe - 2-litre drainage bag with non-return valve and drainage tap | 30 pcs |
| H92011E | TORAMATIC 1-80 | - trio of 80 mm needles with protective cap - 2 one-way valves - 60 ml syringe - 2-litre drainage bag with non-return valve and drainage tap | 30 pcs |
| H92021E | TORAMATIC 1-VER | - 100 mm Veress needle - 2 one-way valves - 60 ml syringe - 2-litre drainage bag with non-return valve and drainage tap | 30 pcs |

KIT FOR PARACENTESIS

| CODE | MODEL | CHARACTERISTICS | PACKAGING |
|---------|------------|--|-----------|
| H92120E | PARASET 2 | - 51 mm packaged cannula needle with 16G diameter - 3-way stopcock - 60 ml syringe - 2-litre drainage bag with non-return valve and drainage tap | 30 pcs |
| H92160E | PARASET 5 | - 51 mm packaged cannula needle with 16G diameter - 3-way stopcock - 60 ml syringe - 5-litre drainage bag with non-return valve and drainage tap | 20 pcs |
| H92200E | PARASET 8 | - 51 mm packaged cannula needle with 16G diameter - 3-way stopcock - 60 ml syringe - 8-litre drainage bag with non-return valve and drainage tap | 20 pcs |
| H92240E | PARASET 12 | - 51 mm packaged cannula needle with 16G diameter - 3-way stopcock - 60 ml syringe - 12-litre drainage bag with non-return valve and drainage tap | 15 pcs |