

COMPACT T SOFT TISSUE BIOPSY


Short notch automatic soft tissue biopsy device

CE
0476

In compliance with MDR
(Regulation EU 2017/745)



- One-handed design simplifies loading and operation whilst allowing for convenient manipulation of the ultrasound transducer
- Programmable: releases the stylet allowing time for tissue to settle in notch (ideal for denser tissue or when there is a desire to verify notch position) and activates the cutting cannula
- Sample notch 10mm for critical procedures
- Thin-wall cannula for larger core
- Colour coded for quick gauge size identification
- Depth markings etched on cannula for precise positioning
- Loaded status indicator
- Plastic safety protection to avoid accidental firing

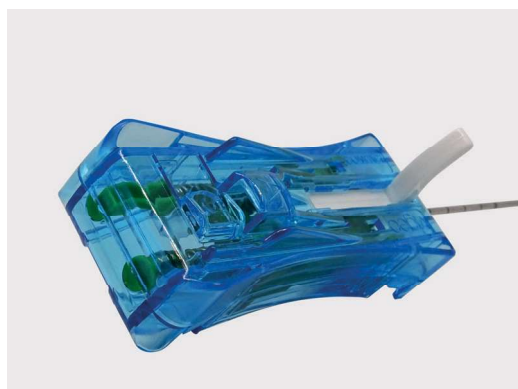

Colour code and plastic safety protection

COMPACT T Qty/Box 10 pcs

CODE	SIZE
● VCT11410	14Gx10cm
● VCT11413	14Gx13cm
● VCT11415	14Gx15cm
● VCT11810	18Gx10cm
● VCT12010	20Gx10cm

On request Vigeo is available to provide needles with different diameter and length