

# PRESSCARE G200M

## PRESSOTHERAPY | 2 OUTPUTS - 2 INFLATION PROGRAMS

MEDICAL DEVICE SUITABLE FOR HEALTH, WELLNESS, BEAUTY AND RELAXATION NEEDS. INDICATED FOR LYMPHEDEMA, PROBLEMS OF VENOUS CIRCULATION, LEGS HEAVINESS, OEDEMAS, SWALLEN LIMBS, LYMPHATIC DRAINAGE, WATER RETENTION, CELLULITE, EMATOMS MESSAGES, VENOUS STASIS. IT CAN BE ALSO USED SIMULTANEOUSLY WITH THREE DIFFERENT ACCESSORIES. IT IS RECOMMENDED TO CONSULT THE DOCTOR BEFORE USING THE DEVICE.

### TECHNICAL FEATURES

<b>Outputs</b>	2
<b>Programs</b>	2
<b>Pressure</b>	0-180 mmHg $\pm$ 10%
<b>Timer</b>	up to 99 minutes
<b>Mains supply</b>	AC 230 V, 50 Hz
<b>Size</b>	185x180x230 mm
<b>Weight</b>	3000 Kg

### EQUIPMENT


- According to the configuration desired
- 1 Power supply unit
- 1 Operating manual

### PURCHASE CONFIGURATIONS

	BOOTS	SLEEVES	ABDOMINAL BAND	SPLITTING CABLE
· PRESSCARE G200M-1	1	/	/	/
· PRESSCARE G200M-2	2	/	/	/
· PRESSCARE G200M-3	2	/	1	1
· PRESSCARE G200M-1B	/	1	/	/



CE  
0598

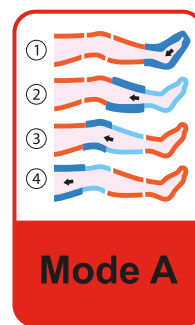

**MODE A:** according to the sequence of the program there will be an inflation and deflation of the boot one sector at a time, from the distal (foot or hand) toward the proximal part.

**MODE B:** according to the sequence of the program there will be an inflation of two chambers at a time, simulating the function of the overlap boots.

# PRESSCARE G200M


## PRESSOTHERAPY | 2 OUTPUTS - 2 INFLATION PROGRAMS

### SIZE ACCESSORIES (4 SECTIONS)


Rev. En\_01\_2022

BOOTS SIZE "S"	(A) LENGTH GROIN - HEEL 72 cm (B) THIGH CIRCUMFERENCE 60 cm with EXTENSION 70 cm (C) ANKLE CIRCUMFERENCE 37 cm with EXTENSION 42 cm
BOOTS SIZE "M"	(A) LENGTH GROIN - HEEL 82 cm (B) THIGH CIRCUMFERENCE cm 67 con ESTENSIONE 77 cm (C) ANKLE CIRCUMFERENCE cm 37 con ESTENSIONE 42 cm
BOOTS SIZE "L"	(A) LENGTH GROIN - HEEL 92 cm (B) THIGH CIRCUMFERENCE cm 74 with EXTENSION 84 cm (C) ANKLE CIRCUMFERENCE cm 40 with EXTENSION 45 cm

SLEEVE	(A) 84 cm - (B) 63 cm (C) SHOULDER CIRCUMFERENCE 43 cm with EXTENSION 53 cm (D) WRIST CIRCUMFERENCE 31 cm with EXTENSION 36 cm
--------	--

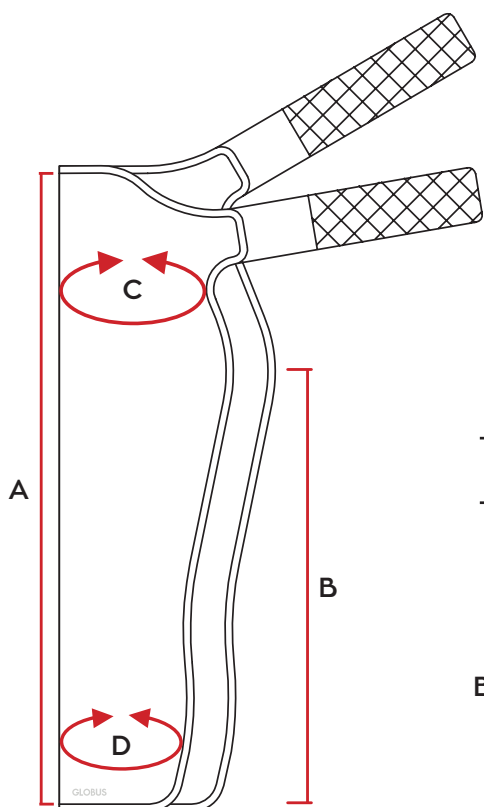

ABDOMINAL BAND	(A) 133 cm - (B) 38 cm ABDOMEN CIRCUMFERENCE from 82 cm to 114 cm
----------------	--

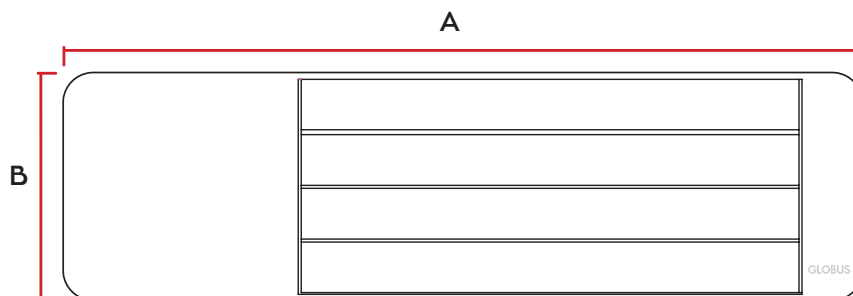

SPLITTER CABLE to use three accessories at the same time      SINGLE CABLE

N.B. All **EXTENSIONS** are to be considered an optional accessory available on request.

The **ACCESSORIES**, depending on the type of configuration, can be optional or included.